MiniSpacer\* Dual Spray Metered Dose Inhaler Adapter Bi-Directional Spray and Cap 22 mm O.D. ISO Connections Accepts plastic and metal pMDI canister stems, including canisters with integrated dose counters.

#### INDICATIONS FOR USE

Model number 1024 is a single patient, disposable device for dispensing pressurized metered dose inhaler (pMDI) medication into a breathing circuit, as prescribed by a physician or other licensed health care practitioner. The device is indicated for patients on a breathing circuit, for whom aerosol medication has been prescribed, in short and long term critical care environments. REF 1024 is intended for use only when connected to a 22 mm fitting.

#### DIRECTIONS FOR USE

**Installing the MINISPACER Adapter in the Patient Circuit**

- Remove the MINISPACER adapter from package.
- Connect the MINISPACER device to the inspiratory limb of the patient circuit, approximately 6 inches back from the patient wye or according to facility guidelines.

CAUTION : Ensure the tubing is secure over MINISPACER connectors when in use.

NOTE : When correctly installed, the MINISPACER MDI port will be vertical and the pMDI canister will point downward during treatment.

- Place the cap firmly on the MINISPACER port.

#### Administering MDI medication via the MINISPACER Adapter

- Uncap the MINISPACER MDI port.
- Insert the MDI canister stem into the MINISPACER port.
- Coordinate MDI actuation with the beginning of inspiration. CAUTION : Support the MINISPACER adapter while depressing the canister.
- Administer the dose as ordered.
- After medication has been delivered as prescribed, remove the pMDI canister and firmly cap the MINISPACER port.
- Monitor the patient response.

NOTE : Discard the MINISPACER adapter when replacing the circuit according to facility guidelines, not to exceed two weeks.

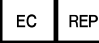
CAUTION: REF 1024 MINISPACER adapter will not increment the counter on pMDI canisters with integrated dose counters surrounding the stem.

WARNING : Do not clean and reuse or sterilize. Potential risks associated with cleaning and reusing, reprocessing and/or sterilization include: cross-contamination or infection resulting from inadequate cleaning, product degradation, and/or failure of the device to perform as intended.

#### Symbol Glossary

	Manufacturer		Consult instructions for use
	Use by Date		Authorized representative in the European Community
	Batch code/lot number		Prescription only. U.S. Federal law restricts this device to sale by or on the order of a physician or properly licensed practitioner.
	CE Mark		Do Not Reuse
	Catalog, reorder or reference number		Caution. The symbol can also be used to mean Attention, see Instructions for Use.

Item No.	Description
 1024	MiniSpacer® Dual Spray MDI Adapter, 22 mm O.D. / O.D.

	Thayer Medical Corporation 4575 S. Palo Verde Rd. Tucson, AZ 85714-1961 USA www.thayermedical.com Made in USA			MDSS GmbH Schiffgraben 41 30175 Hannover, Germany
--	---	---	---	---

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<sup>[4]</sup> trademarks and registered trademarks of their respective companies.

The MINISPACE­R dual-spray adapter was evaluated for particle size and dosing characteristics using an 8 stage cascade impactor. The MINISPACE­R was connected to the cascade impactor with ventilator tubing, wye and an endotracheal (ET) tube as they would be used in a ventilator circuit. Air was drawn through the assembly at a flow rate of 28.3 l/min. The aerosol samples exited the ET tube into the cascade impactor and comprised the Total Delivered Dose.

REF 1024 22 mm O.D. / O.D. – Installed 6 inches upstream of the wye in the inspiratory limb of a 22 mm circuit

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