




# PrimeAire® MDI Holding Chamber

## Features:


- \* 175cc holding chamber
- \* MDI boot adapter accommodates most standard MDI mouthpieces
- \* Standard 22mm OD taper mouthpiece
- \* Dual-Valved system
- \* Low resistance of valve allows patient to breathe normally through the PrimeAire®
- \* Incorporated coaching whistle
- \* Mouthpiece and MDI adapter parts are removable for easy cleaning
- \* Audible coaching whistle

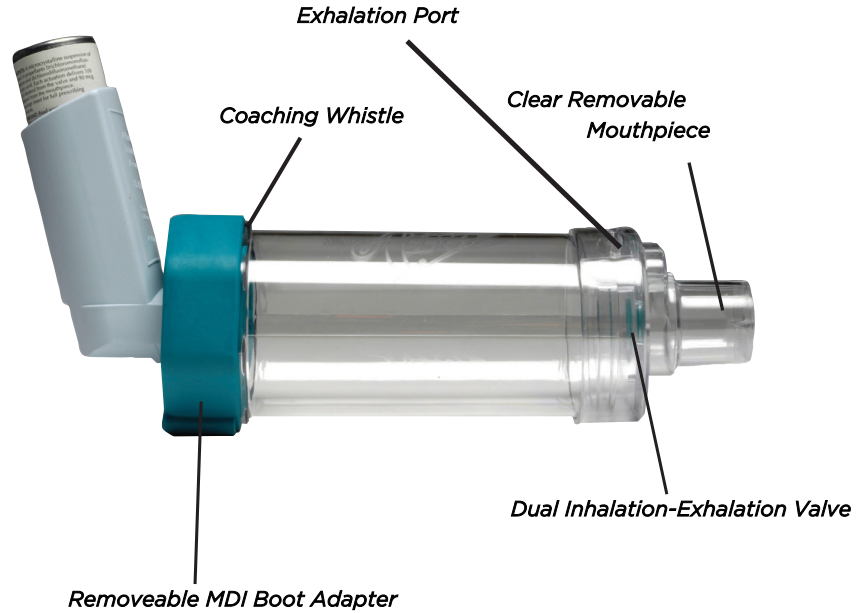

## A clear solution for effective respiratory management

Holding chamber devices, like the PrimeAire®, are commonly used with metered dose inhalers (MDIs) to deliver medication directly to the airways during inhalation. These devices help to improve the direction and deposition of the medication being delivered by an MDI by:

- \* Extending the mouthpiece of the MDI
- \* Retaining the medication in the chamber
- \* Maintaining the particles in suspension

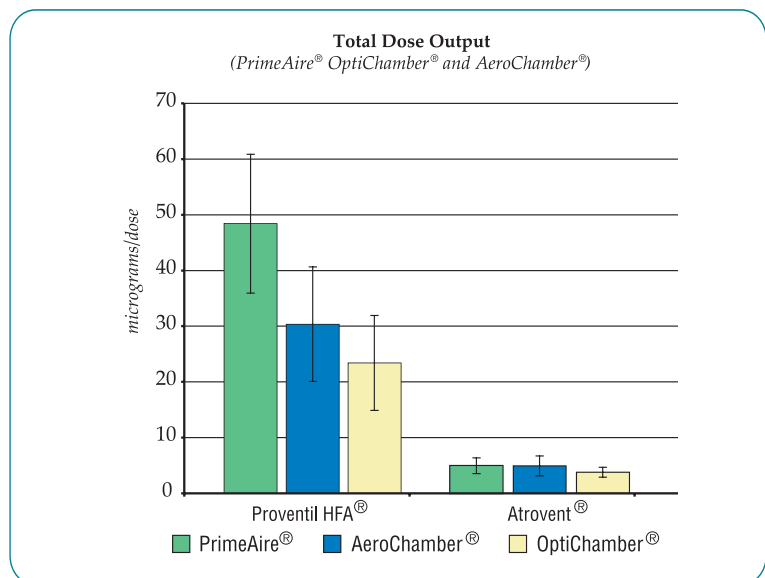

Foss, S., White Paper - Total Dose Output for MDI, Thayer Medical, Tucson, AZ, Fall 1999

4575 S. Palo Verde Rd.  
Ste 337  
Tucson AZ 85714-1961  
www.thayermedical.com

Copyright 2016, Thayer Medical Corporation.  
All rights reserved. PrimeAire® is a registered trademark of Thayer Medical Corporation.  
PUB0302\_rev A.

