



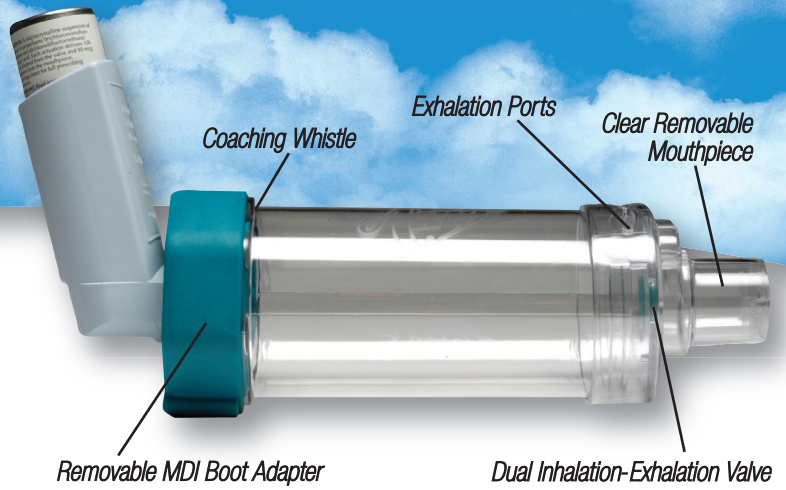

The PrimeAire holding chamber is a clear solution for effective respiratory management.

Holding chamber devices, like the PrimeAire, are commonly used with metered dose inhalers (MDI's) to deliver medication directly to the airways during inhalation. These devices help to improve the direction and deposition of the medication being delivered by an MDI by:

- Extending the mouthpiece of the MDI
- Retaining the medication in the chamber
- Maintaining the particles in suspension

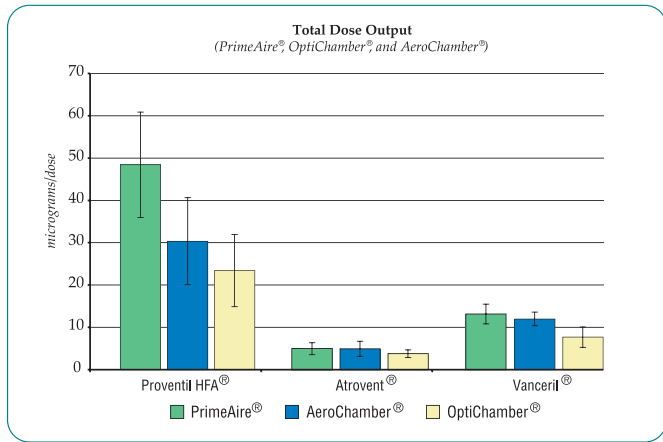

Features:

- 175 cc holding chamber
- MDI boot adapter accommodates all standard MDI mouthpieces
- Compatible with 22mm connection port face masks
- Dual-valve system for safe and efficient inhalation of medication
- Colored dual inhalation-exhalation valve facilitates patient training
- Incorporated coaching whistle alerts the user of improper breathing technique during inhalation

Benefits:

- Boot adapter and mouthpiece easily removed for cleaning
- Reduces common errors associated with MDI use

Cleaning

The PrimeAire is simply cleaned by removing the mouthpiece and MDI boot adapter and washing in warm soapy water. See instructions for use.

Ordering Information

Single Count (Part Number 1114)

For more information, call Thayer Medical at 800-250-3330 or e-mail us at sales@thayermedical.com



THAYER MEDICAL
Thayer Medical
4575 S. Palo Verde Rd., Ste. 337
Tucson, Arizona 85714-1961
(800) 250-3330
www.thayermedical.com

©Copyright 2008, Thayer Medical Corporation. All rights reserved.
PrimeAire® is a registered trademark of Thayer Medical Corporation. PA08-2 rev.A

